



## **User Manual**

**VIA Series Anti-Gravity Treadmill® Product**

This manual covers operation procedures for the following AlterG product:

VIA and VIA X Anti-Gravity Treadmill

Copyright © 2019 by AlterG, Inc. All rights reserved.

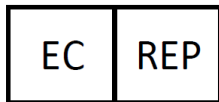
**Note:** The following symbol is used throughout this manual to call attention to Warnings, Cautions or operational procedures that may directly affect the safe operation of the Anti-Gravity Treadmill. Read and understand these instructions and statements before operating the Anti-Gravity Treadmill.



Warning. Cautionary statement or operational procedure that may directly affect the safe operation of the treadmill.

**RoHS**

Compliant product



Obelis s.a  
Bd. General Wahis, 53  
1030 Brussels, Belgium

AlterG and Anti-Gravity Treadmill are registered trademarks in the US by AlterG, Inc.  
Neoprene is a registered trademark of DuPont Corporation.  
Polar is a registered trademark of Polar Electro, Inc.  
Klear Screen is a registered trademark of Meridrew Enterprises.



## Table of Contents

<b>Owner Responsibility</b>	<b>6</b>
<b>AlterG Contact Information</b>	<b>7</b>
<b>Indications and Contraindications for Use</b>	<b>8</b>
Statement of Intended Use	8
Indications for Use	8
Precautions for Use	8
Contraindications for Use	8
<b>Safety: Warnings and Cautions</b>	<b>9</b>
<b>Introduction</b>	<b>11</b>
Consulting a Physician	11
The Importance of Warming Up and Cooling Down	11
<b>Setup and Installation</b>	<b>12</b>
Electrical Requirements	12
Grounding Requirements	12
Location Requirements	12
AlterG Via Anti-Gravity Treadmill® Installations Specifications	12
Transporting Your Anti-Gravity Treadmill	13
Anti-Gravity Treadmill Components	13
<b>Operation</b>	<b>15</b>
Powering On the Anti-Gravity Treadmill	15
Securing the User in the Anti-Gravity Treadmill	15
Starting and Completing a Workout Session	21
Stepping Out of the Anti-Gravity Treadmill	24
<b>Session Metrics</b>	<b>25</b>
Session Metrics	25
Calories	25
<b>Heart Rate Monitoring</b>	<b>26</b>
Heart Rate Monitor	26
<b>Operating HD Video Monitoring System</b>	<b>27</b>
Getting Started	27
Focusing and Positioning Cameras	27
Camera Views	27
<b>Operating Gait</b>	<b>28</b>
Gait Measurements	28
Gait Recording	28
Gait Features	29
Gait Reports	30
<b>Operating Pain</b>	<b>31</b>
Pain Recording	31
Pain Results	31
<b>End-of-Session Reporting</b>	<b>32</b>
Creating End-of-Session Reports	32
Emailing Reports	33
Saving to USB	33

<b>Emergency Extraction Procedures.....</b>	<b>34</b>
Emergency Scenarios.....	34
<b>Facility Settings.....</b>	<b>35</b>
Settings Menu Navigation .....	35
Adjusting Display Layout – GUI Tab.....	35
Manage Emails – GUI Tab.....	36
Tech Support – Support Tab .....	38
Restart App.....	38
Shutdown PC.....	39
Calibrate Deck Weight – Service Tab.....	39
<b>Calibration .....</b>	<b>40</b>
<b>Wi-Fi Connectivity .....</b>	<b>41</b>
<b>Labels, Locations, Interpretation .....</b>	<b>42</b>
<b>Maintenance .....</b>	<b>47</b>
Disinfection .....	47
General Cleaning and Inspection.....	48
Bag Maintenance.....	49
Bag Removal Procedure for Cleaning .....	50
<b>Appendix A: Anti-Gravity Treadmill Specifications.....</b>	<b>54</b>
<b>Appendix B: Options and Accessories .....</b>	<b>56</b>
<b>Appendix C: Troubleshooting .....</b>	<b>57</b>
Repairs.....	57
Air Pressure.....	57
Treadmill .....	57
Leaks.....	57
System Errors .....	57
<b>Appendix D: EMC Statement.....</b>	<b>58</b>
<b>Appendix E: Warranty .....</b>	<b>59</b>

## Owner Responsibility

The AlterG® VIA Anti-Gravity Treadmill® product will perform as described in this manual and by accompanying labels and/or inserts when it is assembled, operated, maintained and repaired in accordance with the instructions provided. The Anti-Gravity Treadmill must be checked periodically as described in this manual. A defective Anti-Gravity Treadmill should not be used. Parts that are broken, missing, plainly worn, distorted or contaminated should be replaced immediately. Should such repair or replacement become necessary, it is recommended that a request for service be made to AlterG, Inc. The Anti-Gravity Treadmill, or any of its parts, should only be repaired in accordance with instructions provided by AlterG, Inc., authorized representatives of AlterG, Inc., or by AlterG, Inc. trained personnel. The Anti-Gravity Treadmill must not be altered without the prior written approval of the AlterG, Inc. Quality Assurance Department.

The owner of this product shall bear the sole responsibility for any malfunction which results from improper use, faulty maintenance, improper repair, damage, or alteration by anyone other than AlterG, Inc. authorized representatives.

The owner/user of the AlterG Anti-Gravity Treadmill has the responsibility to report any serious incident that has occurred in relation to the operation of the device to AlterG and to the nationally recognized competent authority.

Any unauthorized maintenance, repairs or equipment modification activities may void the Anti-Gravity Treadmill Product Warranty.

## AlterG Contact Information

AlterG welcomes your inquiries and comments. If you have any questions or comments, please contact our service and support team.

### **AlterG Headquarters, Service and Support**

48368 Milmont Drive

Fremont, CA 94538 U.S.

+1 510 270-5900

[www.alterg.com](http://www.alterg.com)

# Indications and Contraindications for Use

## Statement of Intended Use

The AlterG VIA Anti-Gravity Treadmill provides unweighting of the user's body weight in addition to normal treadmill functions. The unweighting allows patients and individuals to do standing exercises, walk, or run with reduced impact on their musculoskeletal system.

The VIA is intended for use by individuals with limited mobility who require a start speed of 0.1 mph.

The VIA X is intended for use by individuals with moderate to maximum mobility who are not affected by a start speed of 0.5 mph and a top end speed of 15.0 mph.

## Indications for Use

- Aerobic conditioning
- Sport-specific conditioning programs
- Weight control and reduction
- Gait training and neuromuscular re-education in neurologic patients
- Strengthening and conditioning in geriatric patients
- Rehabilitation following injury or surgery of the lower extremities
- Rehabilitation after total joint replacement

## Precautions for Use

- Cardiovascular disease or respiratory compromise
- Exercise induced asthma or angina
- Acute and chronic back problems
- Ruptured or herniated disc
- Safety and effectiveness in pregnant women have not been established
- Safety and effectiveness for individuals with Functional Independence Measure score of 1 or 2 (dependent of max assist) has not been established
- Treadmill belt does not lock in place. Be sure patient is stable before turning off the treadmill and exiting.
- Heart rate monitor is not accurate.
- Any condition where increased intra-abdominal pressure may be a concern (for example, urinary incontinence, pelvic floor dysfunction, pelvic floor reconstruction, or other conditions)

## Contraindications for Use

- Unstable fracture
- Cardiovascular hypotension
- Deep vein thrombosis



## Safety: Warnings and Cautions

Before using the Anti-Gravity Treadmill, please read this manual. As a physical therapist, trainer, or clinician (the operator of the product), you must understand the safety features and user interface. We want you and your clients or patients (the users of the product) to have a safe and enjoyable exercise experience.



**DANGER: Imminently hazardous situation to be avoided that will result in serious injury or death.**

- Do not modify the Anti-Gravity Treadmill electrical plug. The treadmill comes equipped with a removable power cord with a NEMA 6-20 plug (for 208-240V operation) or a NEMA 5-20 (for 120V operation). You must have the proper receptacle outlet to be compatible. Substitution of this cable to accommodate international receptacles must ensure that the cable is designed to carry 220V at 20A.
- Do not use any electrical adapters. To do so could result in an electrical shock hazard.
- Consult a qualified electrician before using any extension cords. Long extension cords may cause a voltage drop to the Anti-Gravity Treadmill, which may cause it to operate improperly. AlterG provides a 15ft (4.6m) power cord.
- Do not operate the Anti-Gravity Treadmill in wet or damp environments.
- When relevant, do not operate the heart rate monitor transmitter in conjunction with an electrical heart pacemaker or similar device. The transmitter may cause electrical disturbances which can interfere with pacemaker function.
- Always unplug the Anti-Gravity Treadmill before cleaning or servicing.
- Do not soak any part of the Anti-Gravity Treadmill with liquid during cleaning; use a sprayer or damp cloth. Keep all liquids away from electric components. Always unplug the Anti-Gravity Treadmill before cleaning and maintenance.
- Service should be performed by an authorized AlterG technician. Service by non-authorized AlterG technicians will void the warranty. Contact AlterG before you or an electrician attempt any maintenance.
- Do not place any liquids on any part of the Anti-Gravity Treadmill (except in the water bottle holders), including the Anti-Gravity Treadmill running surface.
- Always keep the running surface clean and dry.
- Do not unplug or alter any of the internal wiring on the Anti-Gravity Treadmill after installation.



**WARNING: Potentially hazardous situation to be avoided that could result in serious injury or death.**

- Users must consult with their physicians and obtain a medical exam before beginning any exercise program. This is particularly true if users have any of the following: history of heart disease, high blood pressure, diabetes, chronic respiratory disease, elevated cholesterol, if they smoke cigarettes, are currently inactive, are obese, or have any other chronic disease or physical impairment.
- Users must stop exercising immediately and consult a physician if they feel faint, dizzy, experience chest pains, nausea or any other abnormal symptoms while using the Anti-Gravity Treadmill.



**CAUTION: Potentially hazardous situation to be avoided that may result in minor or moderate injury.**

- Always use the emergency safety lanyard supplied with the Anti-Gravity Treadmill. It should be clipped to the user's article of clothing while exercising. This is an important feature in case the user falls during a workout session.

- Before operating the Anti-Gravity Treadmill, make sure the locking mechanism for the cockpit is engaged and locked all the way into stanchions that frame the Anti-Gravity Treadmill. Failure to do so may allow the cockpit to dislodge during operation, resulting in possible injury to the user.
- Read, understand and test the emergency stop procedure before use.
- Never leave children unsupervised around the Anti-Gravity Treadmill.
- Safety and effectiveness in pregnant women have not been established. Pregnant women or women who may be pregnant should consult their physician before using the Anti-Gravity Treadmill.
- The Anti-Gravity Treadmill must be used under the supervision of a properly trained operator. At no time should a user of the Anti-Gravity Treadmill exercise without appropriate supervision; even if having been previously trained in the proper operation of the device.
- Set up and operate the Anti-Gravity Treadmill on a solid, level surface.
- Do not wear loose or dangling clothing while using the Anti-Gravity Treadmill. Do not store anything (like shorts) inside the Anti-Gravity Treadmill bag.
- Prior to beginning a workout session, check to make sure there is no debris inside the Anti-Gravity Treadmill.
- Keep hands away from the bag and frame structure during inflation to avoid pinching.
- Keep hands away from all moving parts.
- Do not use the VIA Anti-Gravity Treadmill if the user's weight is less than 80lbs (36kg) or greater than 400lbs (180kg). Users who weigh more than 350lbs (159kg) may not be able to unweight to the full 20% of original body weight.
- Care should be taken when the user enters and exits the Anti-Gravity Treadmill. Users should never enter the Anti-Gravity Treadmill while the treadmill surface is moving. Ensure that the emergency stop safety magnet is attached to the treadmill so that the treadmill belt is locked and will not move when the user steps on the surface. Make sure the user holds onto the cockpit or handrails whenever practical to support their body.
- Make sure that the user is fully zipped into the bag before beginning the workout session and that the cockpit is adjusted at the correct height and locked in place.
- Make sure the user wears proper athletic shoes, such as those with rubber or high-traction soles. Do not allow shoes with heels or leather soles. Make sure no stones or sharp objects are embedded in the soles of the shoes.
- As with any treadmill workout, make sure the user includes a cool-down phase at the end of the user's workout session. Make sure they return to full body weight and exercise moderately before stopping. Avoid abruptly ending or pausing the workout session while the user is at reduced body weight or at high speed.
- The safety and integrity of the Anti-Gravity Treadmill can only be maintained when the Anti-Gravity Treadmill is regularly examined for damage and wear and is properly repaired. It is the sole responsibility of the user/owner or facility operator to ensure that regular maintenance is performed. Worn or damaged components must be replaced immediately, and the Anti-Gravity Treadmill removed from service until the repair is made. Only manufacturer supplied or approved components should be used to maintain and repair the Anti-Gravity Treadmill.



**CAUTION: AlterG uses Neodymium magnets to hold the bag in place near the user. Neodymium magnets are a member of the Rare Earth magnet family and are the most powerful permanent magnets in the world. The VIA product uses up to four N42 Neodymium magnets on each side of the bag. Each of these magnets can have as much as 22 lbs. (10kg) of pulling force and a magnetic field of 13,200 Gauss. Fortunately, the magnetic field drops very quickly from the magnet such that, at a distance of 3 inches (7.6cm), the magnetic field is less than 100 Gauss in total. Caution is advised if a user has a heart pacemaker, insulin pump, or other electromagnetic products in use during a session within the AlterG Anti-Gravity Treadmill. Always check with your healthcare professional prior to use of this product.**

# Introduction

## Consulting a Physician

Anyone considering an exercise program or an increase in activity should consult a physician. It is highly recommended that users follow the guidance of their physician before and during an exercise program or any other increase in physical activity if they:

- Have heart disease, high blood pressure, diabetes, chronic respiratory disease, or elevated cholesterol
- Smoke cigarettes
- Are currently inactive, are obese, or have any other chronic disease or physical impairment, or if there is a history of such disease in their family



### **Neodymium magnets can affect pacemakers.**

The strong magnetic fields near a neodymium magnet can affect pacemakers, ICDs and other implanted medical devices. Many of these devices are made with a feature that deactivates it with a magnetic field. Therefore, care must be taken to avoid inadvertently deactivating such devices.

## The Importance of Warming Up and Cooling Down

It is important that users gradually warm up, cool down, and incorporate a series of stretches prior to and at the end of each work out. Stretching encourages the necessary flexibility to help prevent sore muscles and injury during daily activities.

Do not abruptly end the workout session on the VIA Anti-Gravity Treadmill. The user's full body weight should be restored slowly and should include a few minutes of walking at full body weight and low intensity before stopping the workout session.

## Setup and Installation

The VIA Anti-Gravity Treadmill will be installed by an AlterG qualified technician after delivery. Please make sure that you inspect the Anti-Gravity Treadmill upon delivery for any damage that may have occurred during transportation. Take photographs and report any damage immediately to the shipping company and AlterG. When you sign for the shipment of your Anti-Gravity Treadmill, you are taking responsibility for any damage that may occur before installation.

## Electrical Requirements

The recommended electrical power connection for the VIA Anti-Gravity Treadmill is a dedicated 20 ampere, 220VAC @ 50/60 Hz circuit with ground. The plug supplied with the Anti-Gravity Treadmill is designated by the NEMA configuration system as 6-20P. The corresponding receptacle for the plug is a NEMA 6-20R.

An appropriate plug for the country and facility in question should be wired as follows:

Blue Conductor: Neutral  
Brown Conductor: Line  
Green/Yellow Conductor: Ground

## Grounding Requirements



**WARNING:** To avoid the risk of electric shock, this equipment must only be connected to a supply main with protective earth.

The VIA Anti-Gravity Treadmill must be grounded electrically. If there is an electrical malfunction, grounding provides a path of least resistance for electric current to reduce the risk of electric shock. The Anti-Gravity Treadmill power cord includes a grounded plug. This plug must be plugged into an appropriate receptacle -that is properly installed and grounded in accordance with the current National Electrical Code as well as all local codes and ordinances. If you are at all unsure of these requirements, contact AlterG or a qualified electrician.

## Location Requirements

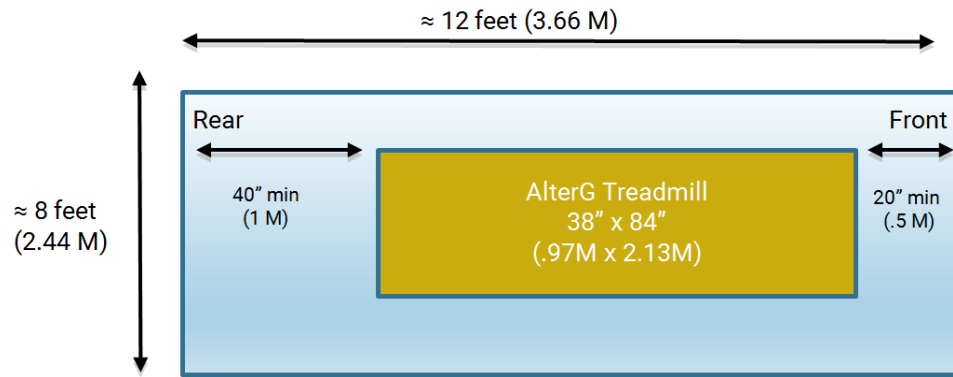
The VIA Anti-Gravity Treadmill measures 85" (217cm) × 38" (97cm) and weighs almost 550lbs (249kg). It needs to be placed on a structurally sound surface. If it is used above ground level, it is advisable to place it near the corner of the room or where the floor will be strongest to ensure maximum support during high-speed, high-impact use. The surface should be level to ensure minimal flexing of the Anti-Gravity Treadmill frame. Do not place the Anti-Gravity Treadmill directly on thick carpeting because it may interfere with air-valves located on the underside of the treadmill. The Anti-Gravity Treadmill needs to be within 12ft (3.7m) (from the front of the treadmill) of the proper electrical outlet. Check with a qualified electrician or AlterG if you plan on extending the cord in any way. Make sure you leave at least 24" (61cm) on either side of the treadmill to allow the bag to expand during inflation. Allow at least 40" (1m) behind the treadmill to accommodate a user getting in and out safely, and at least 20" (.5m) in front of the unit to allow access to electronics and power on/off switch.

We recommend an area at least 12ft (3.7m) long by 8ft (2.4m) wide to provide adequate space for operation and user access. Also, check ceiling height to ensure that users will not hit their heads on the ceiling while running. The Anti-Gravity Treadmill surface is 6" (15cm) off the floor. An 8ft (2.4m) ceiling may be too low for taller users.

## AlterG VIA Anti-Gravity Treadmill® Installations Specifications

- Maximum height of treadmill: 74" (1.88M)
- Width of treadmill: 38" (0.97M)

- Length of treadmill: 85" (2.17M)
- Ceiling height minimum: 96" (2.44M)
- Product Weight is ~550lbs (249Kg)
- **Important note: Place the front within 8 feet (2.44M) of the power outlet and leave at least 20 inches (0.5M) clearance in the front of the unit for access to electronics and on/off switch.**

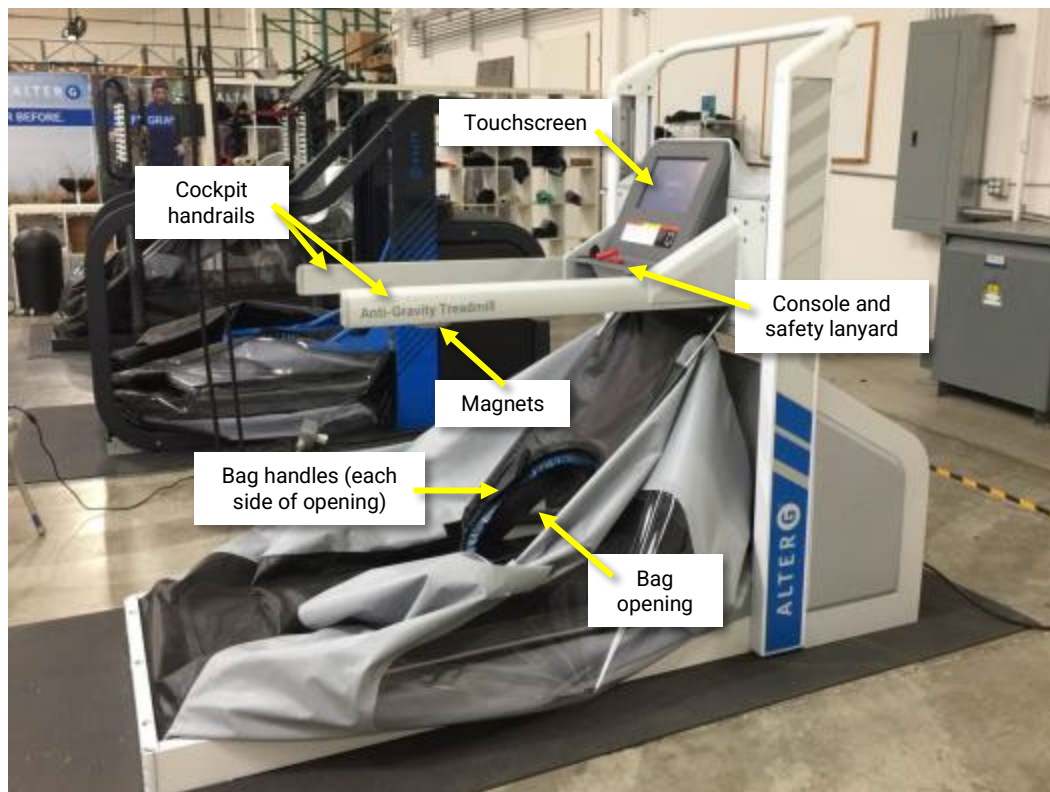


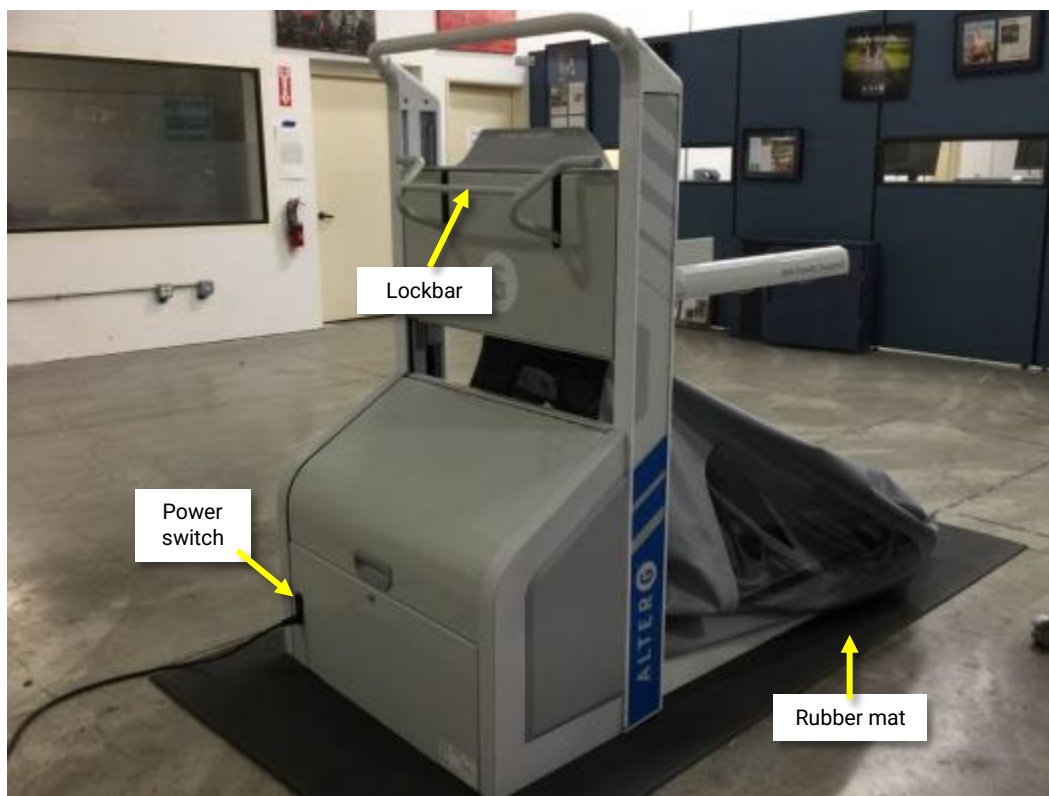
## Transporting Your Anti-Gravity Treadmill

Contact AlterG if you plan on moving your Anti-Gravity Treadmill. Damage sustained by improperly moving the Anti-Gravity Treadmill will not be covered by your warranty.

## Anti-Gravity Treadmill Components

Key components for operating the Anti-Gravity Treadmill are labeled in the photos below.





## Operation

**Note:** The following operating instructions are performed by, or at the direction of, a clinician.

### Powering On the Anti-Gravity Treadmill

1. The power switch is located at the front of the VIA Anti-Gravity Treadmill. Switch it on to start the system and to boot up the software. As the software loads, avoid touching the screen because it may interfere with the boot-up process.

**Note:** Do not allow the user to stand on the Anti-Gravity Treadmill belt during boot-up. If the user inadvertently puts weight on the belt during boot-up, an error message is displayed. Have the user stand outside the unit until the boot sequence is complete.



2. The system is ready when the Welcome screen appears.

**Note:** If you have not read this manual, it is recommended that you QUIT now or seek appropriate instruction from a qualified operator. By tapping **START SESSION**, you acknowledge that you have read and understand this manual.



### Securing the User in the Anti-Gravity Treadmill

The procedures in this section describe:

- Putting on the shorts
- Stepping into the Anti-Gravity Treadmill
- Zipping into the bag
- Setting the cockpit height
- Attaching the safety lanyard

#### Putting on the Shorts

Your Anti-Gravity Treadmill comes with customized neoprene compression shorts that ensure an airtight seal between the user's body and the bag opening of the treadmill. It is recommended that the user put the shorts on before stepping into the treadmill.

Have the user select a size that is snug but not uncomfortable and make certain the tag is at the user's back and on the inside of the shorts. For a comfortable fit and to prevent bunching, it is recommended that the user wear a pair of running shorts or tights under the shorts.

**Note:** The shorts are latex-free (90% Neoprene, 10% Urethane coated Nylon).





## Stepping into the Anti-Gravity Treadmill



**CAUTION:** Before the user steps into the Anti-Gravity Treadmill running surface, ensure that the safety magnet is in place on the console. If the safety magnet is not in place, the treadmill assumes there is a safety problem and disengages the running surface.

Before entering the treadmill, ensure that the cockpit is at the lowest level. Entering with the treadmill with the cockpit in a higher position may cause the cockpit to move, potentially causing the user to lose balance.

Enter the treadmill from the back and step into the bag opening.

It is fine to step on the bag, but make sure that no rocks or sharp objects are embedded in the soles of shoes that could mar or damage the bag.

**Note:** The bag fabric does not contain latex.



## Setting the Cockpit Height

1. Face forward and pull the lock bar toward you to disengage the cockpit.





---

A “Cockpit Lock” message appears on the screen to indicate that the cockpit is unlocked.



2. Grasp the handrails and slide the cockpit up or down to set it at a height that places the zipper at the iliac crest (the top, outer edge of the pelvic bone felt just below waist level). For additional trunk support, the cockpit can be set higher.

**Note:** Make sure there is enough room above the handrails to allow your arms to swing comfortably during the workout session.



**Note:** Inseam heights are posted on the side stanchions. These can be used to mark the cockpit height for different users.

For more information, see the description of the inseam height labels in “Labels, Locations, Interpretation”.



3. With the cockpit at the correct height, make sure it “clicks” into place, and then push the lock bar forward (away from you) until it is fully engaged.



---

The “Cockpit Lock” message disappears, and the Welcome screen appears.



- 
4. Grasp the handle on each side of the bag and lift the bag straight up to the handrails.



- 
5. Pull up until both sides of the bag snap to the magnets on the underside of the handrails.





**CAUTION:** Never operate the Anti-Gravity Treadmill without the cockpit engaged in a locked position. This is an essential safety item that prevents the cockpit from inadvertently disengaging during an exercise session.

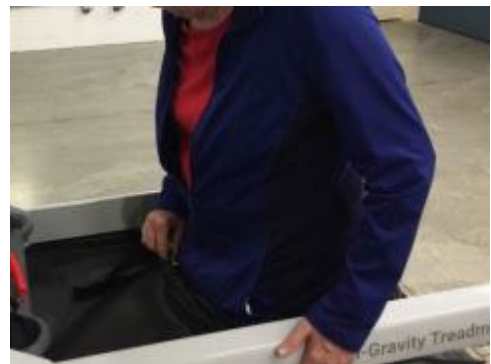
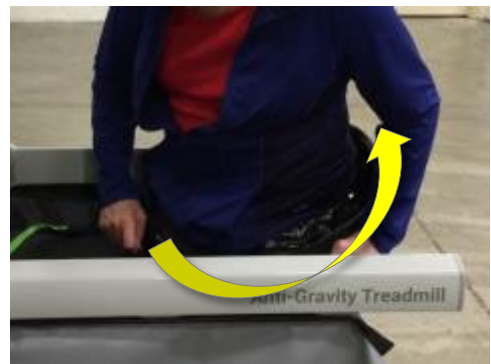
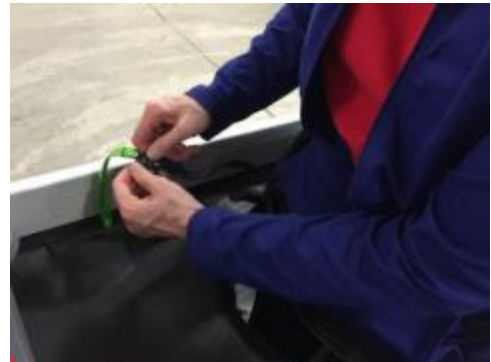
Never attempt to move the cockpit while the bag is inflating or when it is fully inflated. If the bag needs to be repositioned while the user is exercising, stop the session. You can then change the height of the cockpit.

---

## Zippering into the Bag




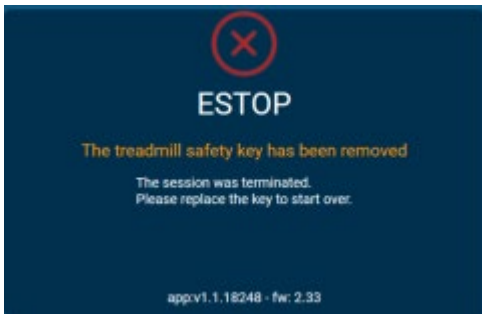
When the cockpit is in place, zip into the Anti-Gravity Treadmill. The zipper is started at the front and center of the body and zipped counter-clockwise all the way around until it returns to overlap in the front. Make sure that the zipper is completely closed to complete the seal.

The zipper provides a means of quickly attaching and detaching the user from the bag, simplifies hygiene, and provides a custom fit for users of all sizes.



## Attaching the Safety Lanyard

Always use the magnetic safety lanyard supplied with the Anti-Gravity Treadmill. If you lose the safety magnet, order a replacement from AlterG. As a precaution, the treadmill will not operate without the safety magnet directly over the Emergency Stop label.

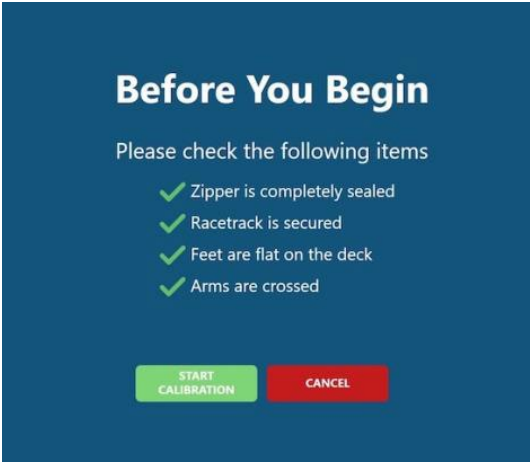
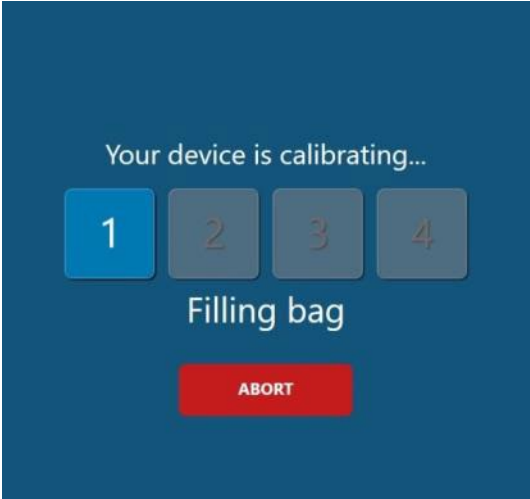
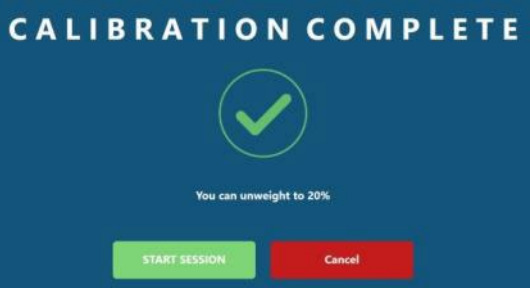
Operator or User Action	Screen
1. Slip the wrist strap of the lanyard around one wrist.	
2. Place the lanyard safety magnet directly over the circular indicator on the console labeled "Emergency Stop"; otherwise, the Anti-Gravity Treadmill will not operate.   <b>CAUTION:</b> Never attempt to defeat this critical safety feature by fastening the lanyard to the bag, structure of the Anti-Gravity Treadmill, or anywhere else other than to the user's wrist.  Should the user fall while exercising, the magnet is pulled off the console and power to the treadmill is cut; air pressure is released; the treadmill running surface disengages from its drive and can move freely. The ESTOP message is displayed on the screen.	 
3. Continue to "Starting and Completing a Workout Session".	

## Starting and Completing a Workout Session

The procedures in this section describe:

- Starting a workout session
- Using the touchscreen features and controls
- Adjusting Anti-Gravity Treadmill speed

### Starting a Workout Session

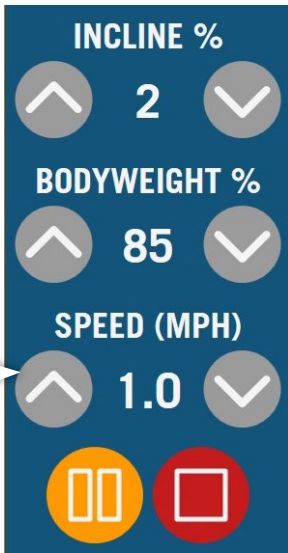
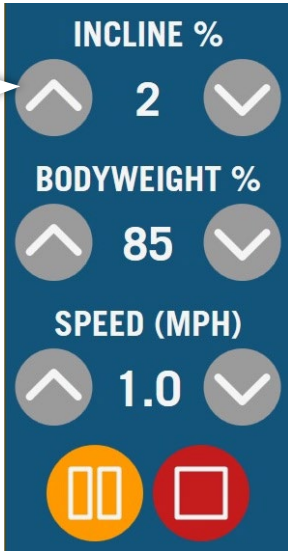
Operator or User Action	Screen
<p>1. PREPARATION FOR EXERCISE.</p> <p>Before you begin, stand still on the surface of the treadmill belt. Do not hold onto or support yourself on any part of the system structure. The system weighs you prior to exercise and the treadmill must support your full body weight.</p>	 <p>The screen has a dark blue background. At the top, the text 'Before You Begin' is in white. Below it, 'Please check the following items' is in a smaller white font. There are four green checkmarks, each followed by a white text item: 'Zipper is completely sealed', 'Racetrack is secured', 'Feet are flat on the deck', and 'Arms are crossed'. At the bottom, there are two buttons: a green 'START CALIBRATION' button and a red 'CANCEL' button.</p>
<p>2. BEGINNING THE EXERCISE SESSION.</p> <p>Tap <b>START CALIBRATION</b> to begin the exercise session. The Anti-Gravity Treadmill will proceed with a calibration routine that allows the system to determine the relationship between enclosure pressure and your body weight. You will feel the pressure in the bag change as the system calibrates. It is best to cross your arms while the routine runs to avoid touching any part of the structure and ruining the calibration.</p>	 <p>The screen has a dark blue background. At the top, the text 'Your device is calibrating...' is in white. Below it are four numbered boxes: '1' (blue), '2' (grey), '3' (grey), and '4' (grey). Below the boxes, the text 'Filling bag' is in white. At the bottom, there is a red 'ABORT' button.</p>
<p>You will be promptly notified when the calibration sequence is complete, tap <b>START SESSION</b> to start your workout session. Inactivity will cancel your completed calibration.</p> <p>This screen also displays your lower body weight limit.</p>	 <p>The screen has a dark blue background. At the top, the text 'CALIBRATION COMPLETE' is in white. Below it is a green checkmark inside a circle. Below the checkmark, the text 'You can unweight to 20%' is in a smaller white font. At the bottom, there are two buttons: a green 'START SESSION' button and a red 'Cancel' button.</p>

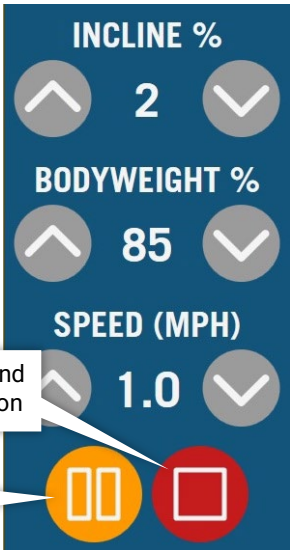
The controls on the Main Screen are shown below.



Operator or User Action	Screen
<p>3. ADJUST YOUR BODY WEIGHT.</p> <p>Adjust your body weight with the up arrow and down arrow button controls. Tapping up arrow increases your body weight and tapping down arrow decreases your body weight. Your body weight is displayed as a percentage of your full body weight. You can adjust your body weight at any time during your exercise session.</p> <p><b>CAUTION:</b> At 40% body weight and lower, you can become unstable if you jump or perform any other activity besides walking or running. Reduce your body weight percentage slowly so you can become accustomed to the new sensation and adjust your gait mechanics accordingly.</p>	<p>The close-up shows the BODYWEIGHT % control with an up arrow and a down arrow. The current value is 85%. Callouts indicate that tapping the up arrow increases body weight and tapping the down arrow decreases body weight. A note below the control states: '% of your full body weight'.</p>



Operator or User Action	Screen
<p>4. SELECT TREADMILL DIRECTION.</p> <p>When the treadmill turns on, it is programmed for forward ambulation. To walk backwards, tap down arrow until it displays negative numbers.</p>	
<p>5. ADJUST VIA TREADMILL SPEED.</p> <p>Adjust treadmill speed using the up arrow and down arrow button controls. In the forward direction, speed will start at 0.1 mph (0.16km/hr) and increase in 0.1mph (0.16km/hr) increments for each button tap up. When the treadmill is programmed in reverse, speed will start at 0.1 mph (0.16km/hr) and increase in 0.1mph (0.16km/hr) increments.</p> <p>Note, recent software changes now allow for full 1 MPH (1.6Km/Hr) increments speed changes. This feature can be turned on/off through set-up in the facilities screen when using software version 1.1.21084 or higher.</p>	
<p>6. ADJUST TREADMILL INCLINE.</p> <p>Tap up arrow to increase incline and tap down arrow to decrease incline. Holding the button down for more than 2 seconds will cause incline to change at a more rapid rate.</p> <p><b>Note:</b> Treadmill incline cannot be adjusted greater than 5% if speed is set for less than 0.5mph (0.8km/hr).</p>	

Operator or User Action	Screen
<p>7. ENDING OR PAUSING THE EXERCISE SESSION.</p> <p>Before stopping your session, return the treadmill to 0% incline. Tap <b>STOP</b> to end the exercise session.</p> <p>To place the treadmill in pause mode, tap <b>PAUSE</b>. The tread belt will stop, but all statistical information will be preserved. Body weight support is maintained while the treadmill is paused. To resume the exercise session, tap <b>PAUSE</b>.</p>	 <p>Tap STOP to end exercise session</p> <p>Tap PAUSE to pause treadmill</p>

## Stepping Out of the Anti-Gravity Treadmill

1. Stop the session and wait for the treadmill to come to a complete stop.
2. Wait for the bag to completely deflate.
3. The user can then:
  - Remove the safety wrist strap.
  - Unzip the shorts from the bag.
  - Disengage the cockpit locking mechanism by pulling the lockbar toward them.
  - Gently lower the cockpit.
  - Exit the treadmill. Turn around, step out of the bag opening, and carefully step off the back of the treadmill.





## Session Metrics

The session screen gives you a graphical overview of body weight, speed, and inclines. It also includes information on your average speed, calories burned, and total distance.

## Session Metrics

The session metrics, by default, is displayed on the top portion of the screen. The session time, calories burned, pace, distance traveled, and heart rate are displayed. Alternative layouts can be set in the facility setting. For more information, see “Adjusting Display Layout – GUI Tab”.



## Calories

The calorie count displayed considers your body weight percentage and treadmill speed, incline, and distance. As a result, the calories burned reading is more accurate when compared to exercise equipment that does not include these parameters.

Reference <http://42.195km.net/e/treadsim/> for more information on how the VIA calculates calories burned.

# Heart Rate Monitoring

## Heart Rate Monitor

The AlterG Treadmill display is designed to receive a user's heart rate in conjunction with the use of a Polar® (Chest Strap) Heart Rate Monitor.<sup>1</sup> Polar chest straps can be purchased at most popular sporting goods stores or online. In order for the screen to correctly display a user's heart rate, the receiver within the display must obtain a stable heart rate signal from the Polar transmitter. The Polar Heart Rate System consists of two main elements: 1) the sensor/transmitter and 2) the receiver within the AlterG Treadmill display.

## How to Wear the Chest Strap Transmitter

The Sensor/Transmitter is worn just below the chest and at the top of the abdomen, directly on bare skin (not over clothing). The transmitter should be centered below the pectoral muscles. Once the strap is secured, pull it away from the chest by stretching the band, and moisten the conductive electrode strips with plain water. The transmitter operates automatically while you are wearing it; it does not operate while it is disconnected from your body. However, as moisture may activate the transmitter and salt buildup from sweat can be a problem, rinse the transmitter with water and wipe it dry after use. The chest band is washable. After you have detached the transmitter, wash the band in warm water using mild soap and rinse thoroughly in clean water.

## The Receiver

You must be within two and a half feet of the receiver for the signal to be received. Please take note that your transmitter may fluctuate erratically if you are too close to other Polar equipment. Maintain at least a three-foot distance between other Polar units.

**NOTE: Erratic heart rate reception may occur if the Polar Monitor is too close in proximity to strong sources of electromagnetic radiation, such as television sets, Personal Computers, electric motors and some other types of fitness equipment. Only one transmitter should be used inside the range of any one receiver as the receiver may pick up several signals simultaneously, causing an inaccurate readout.**




**CAUTION: The heart rate obtained while on the AlterG Treadmill is for reference only. No medical treatment or diagnosis should be based on the Polar Monitor.**

---

<sup>1</sup> Polar® is a registered trademark of Polar Electro, Inc. Lake Success, NY

# Operating HD Video Monitoring System

## Getting Started

Operator or User Action	Screen
To view the camera, tap <b>VIDEO</b> on the navigation bar.	

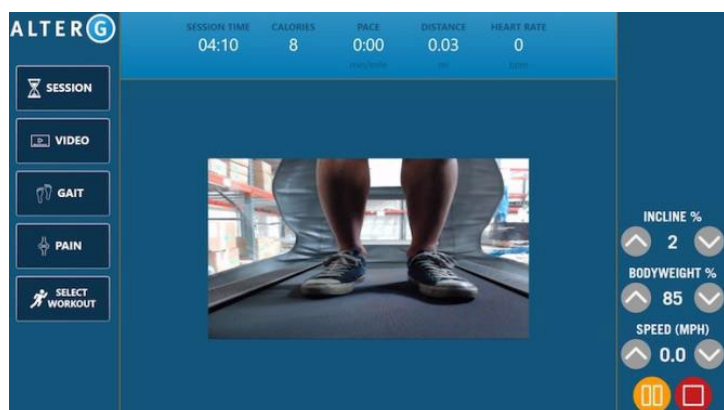
## Focusing and Positioning Cameras

There is a single camera on the VIA Anti-Gravity Treadmill. The front camera is mounted on the inside of the bag and focused at the factory before shipment.



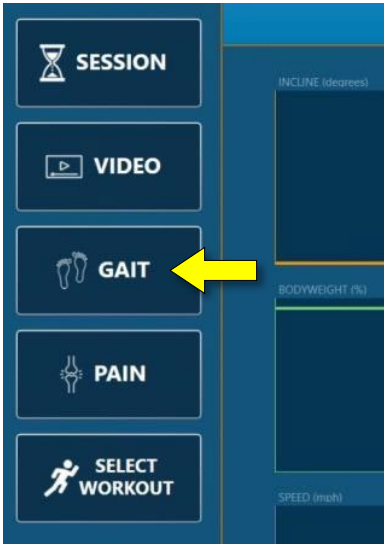
## Camera Views

The front camera is permanently affixed inside the treadmill bag. The video feed will be displayed on the touch monitor on the main display window.




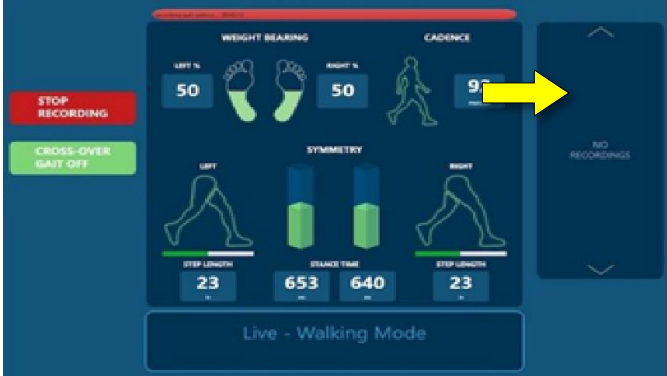
# Operating Gait

## Gait Measurements

Operator or User Action	Screen
<p>The Anti-Gravity Treadmill measures the following:</p> <ul style="list-style-type: none"><li>• Weight bearing symmetry</li><li>• Cadence</li><li>• Stance time symmetry</li><li>• Step length symmetry</li></ul> <p>To view Gait, tap <b>GAIT</b> on the navigation bar.</p> <p>Gait measurements are a display of the last 3 steps on each foot.</p> <p>NOTE: Gait features are intended for use at moderate walking speeds only.</p>	 A screenshot of the treadmill's control interface. On the left is a vertical navigation bar with five buttons: 'SESSION' (hourglass icon), 'VIDEO' (play icon), 'GAIT' (foot icon, highlighted with a yellow arrow), 'PAIN' (joint icon), and 'SELECT WORKOUT' (runner icon). On the right, a main display area shows 'INCLINE (degrees)' at the top, followed by 'BODYWEIGHT (%)' with a green progress bar, and 'SPEED (mph)' at the bottom.


## Gait Recording

Operator or User Action	Screen
<p>Users can record as many gait reports as they want. There is no time limit to the length of the recording. To record a gait report, navigate to the gait page and tap <b>RECORD GAIT DATA</b>.</p> <p>A red bar will appear at the top of the screen signaling that gait is currently being recorded. Tap <b>STOP RECORDING</b> when you want to stop recording.</p> <p>You can adjust any of the treadmill controls while recording (body weight %, incline, and speed). We recommend recording different reports so that you can see the effect of body weight, incline, or speed on gait symmetry.</p>	 A screenshot of the 'Gait Recording' screen. On the left, there are two green buttons: 'RECORD GAIT DATA' (highlighted with a yellow arrow) and 'CROSS-OVER GAIT OFF'. The main area displays several metrics: 'WEIGHT BEARING' with 'LEFT %' at 33 and 'RIGHT %' at 67; 'CADENCE' at 126; 'SYMMETRY' with 'LEFT' step length at 6, 'STANCE TIME' at 320, 'RIGHT' step length at 14, and '790' in the center. At the bottom, it says 'Live - Walking Mode'. On the right, a vertical bar shows 'NO RECORDINGS'.


Operator or User Action	Screen
When you have completed your recording, reports are automatically listed to the right of the data screen.	

## Gait Features

There is one additional feature for different use cases on the Gait Screen.

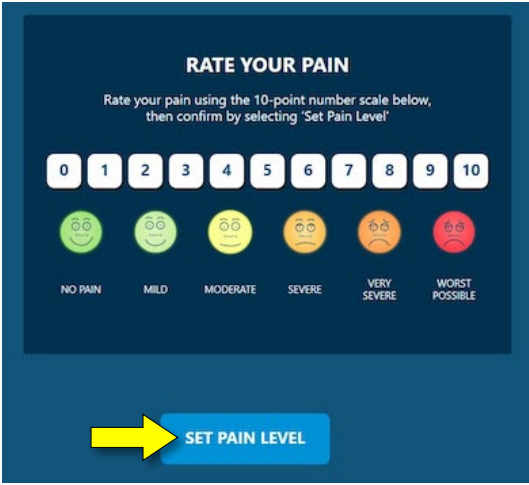
Operator or User Action	Screen
<p><b>Cross-over gait.</b> If you have a patient with scissoring gait (foot crosses the midline when planted), tap <b>CROSS-OVER GAIT</b>. The button is highlighted in orange when in effect. This allows the software to identify the right foot when the left load cells detect weight, and the left foot when the right load cells detect weight.</p>	

## Gait Reports


Operator or User Action	Screen
<p>To view previous Gait Reports, navigate to the Gait screen. All reports are listed to the right of the Gait Data.</p> <p>Select which gait report(s) you want to review from the thumbnails.</p> <p>A summary of that data session will display on the screen. This is an average of each of the gait measurements over the duration of your recording. (that is., if you recorded for 30 seconds, this is a summary of your gait symmetry over those 30 seconds).</p> <p>The thumbnails are labeled with the starting body weight % and session time of your recording to allow you to easily locate the report you want to view.</p> <p><b>Note:</b> You will be able to email or save gait reports at the end of your session (for more information, see "Creating End-of-Session Reports").</p>	 <p>The screenshot shows the Gait screen interface. On the left, there are two green buttons: "RECORD GAIT DATA" and "CROSS-OVER GAIT OFF". The main area displays several metrics: "WEIGHT BEARING" with "LEFT % 51" and "RIGHT % 49", "CADENCE" with "150 steps/min", "SYMMETRY" with "LEFT 35", "DISTANCE TIME" with "407" and "393", and "RIGHT 34". At the bottom, it says "Live - Running Mode". On the right side, there is a vertical list of reports, each with a thumbnail and a yellow arrow pointing to it. The top report shows "85% BODYWEIGHT 00:00".</p>

# Operating Pain

## Pain Recording

Operator or User Action	Screen
<p>To record pain levels during the exercise session, tap <b>PAIN</b> on the navigation bar.</p> <p>Pain can be recorded on a scale of 0-10 (0 = no pain, 10 = maximum pain). To record pain, select a number from 0-10 and tap <b>SET PAIN LEVEL</b>.</p> <p>There are no prompts during the session to remind you to input pain. It is up to the therapist to ask the patient to set their pain levels throughout the session.</p> <p>After you have set the pain level, the View Results page is displayed.</p>	

## Pain Results

Operator or User Action	Screen
<p>To review previous pain recordings, tap <b>VIEW RESULTS</b>.</p> <p>A graph displays all recorded pain results. The graph displays pain in relation to approximate body weight %. (for example, If you recorded a pain level of 8 at 93% body weight, the number 8 is listed between 90% and 100%).</p> <p>Tolerance is also displayed. This is a measurement of the body weight % where the patient felt the least amount of pain.</p> <p><b>Note:</b> Pain results are automatically included when you email or save your end-of-session reports (for more information, see “Creating End-of-Session Reports”).</p>	

## End-of-Session Reporting

The end-of-session report is a summary of your session details. At the end of your exercise session, you can email or save a copy of your session report, gait data, and pain measurements captured during your session. It will be sent to you by email or saved to USB in PDF format. There will be no patient or athlete name on the report.

Any session reports not emailed or saved at the end of the session will be deleted.

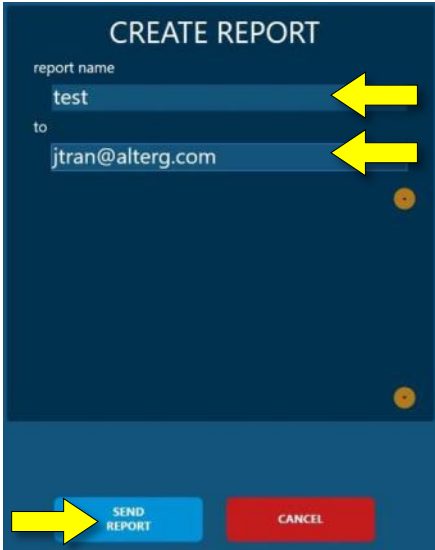
Estimated wait time to email reports is approximately 20 seconds, depending on your Wi-Fi connection. Please be patient and confirm all data has been sent before exiting the treadmill.

## Creating End-of-Session Reports


Operator or User Action	Screen
1. To create a report, tap <b>CREATE REPORT</b> .	
2. Select the reports you want to email or save by tapping the thumbnails you want. Selected thumbnails are highlighted. You can send as many reports as you want.	
3. After you have selected all of your items, tap <b>SEND BY EMAIL</b> , or tap <b>EXPORT TO USB</b> to save to USB.	



## Emailing Reports

Operator or User Action	Screen
<ol style="list-style-type: none"><li>1. If you select Send by Email, enter a report name in the Create Report screen. We recommend a non-identifiable patient ID (for example, 12345), or the session time and date (for example, 181205 02:00). This will help you identify the patient. We recommend emailing the information to yourself and later attaching it to the patient's EMR.</li><li>2. Enter your email address and tap <b>SEND REPORT</b>.</li></ol>	

## Saving to USB

Operator or User Action	Screen
<ol style="list-style-type: none"><li>1. If you select Export to USB, enter a report name in the Create Report screen. We recommend a non-identifiable patient ID (for example, 12345), or the session time and date (for example, 181205 02:00). This will help you identify the patient.</li><li>2. Tap <b>SAVE REPORT</b>.</li></ol> <p><b>Note:</b> You can save as many reports as you want, if there is enough space on your thumb drive.</p>	

# Emergency Extraction Procedures

In the unlikely event that a patient becomes immobilized in the Anti-Gravity Treadmill, please refer to the urgent steps below to secure the situation before calling +1510.270.5900 or referring to [alterg.com](http://alterg.com).

## IMPORTANT NOTES:

1. Do not unlock the cockpit until you read further.
2. The treadmill belt on the VIA series does not lock in place, so secure it by holding your foot against it at the rear of the treadmill.

## Emergency Scenarios

There are two likely scenarios in which a patient may require urgent help with getting out of the Anti-Gravity Treadmill:

### Scenario 1

The patient falls or becomes fatigued to the point of inability to exit the treadmill before the air chamber is inflated with air.

- Option 1: If the patient can stand, stop the treadmill and secure the belt with your foot. Help the patient exit the treadmill.
- Option 2: Prevent the belt from moving by securing it with your foot. Remove the bag nuts at the treadmill base.
- Option 3: Prevent the belt from moving by securing it with your foot. Cut the air chamber and help the patient exit the treadmill.


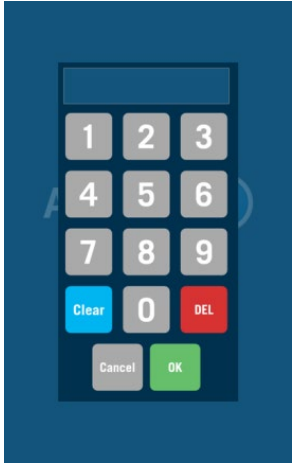
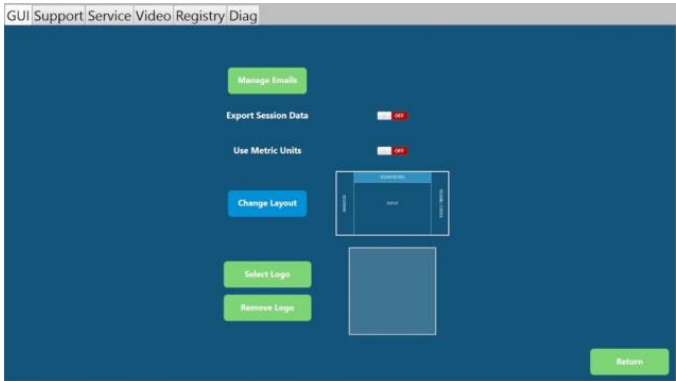


### Scenario 2

The patient becomes fatigued to the point of inability to exit the treadmill after the air chamber is inflated.

Reduce the body weight percentage to less than 40% and have the patient sit down in place until they regain their strength and can exit on their own.

## Facility Settings

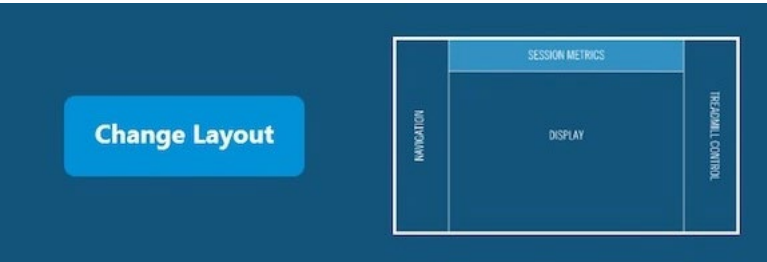
Operator or User Action	Screen
1. To access all setting screens, return to the home screen and double-tap on the upper-left corner of the screen.	
2. A keypad will pop up. Enter code <b>5900</b> .	
Upon entering the settings menu, you will be greeted with a few configurations and options.	

## Settings Menu Navigation

Your facility will have limited access to the Settings navigation bar to aid in troubleshooting and diagnostics. You will be granted access to the GUI, Support, and Service and other highlighted tabs. Inaccessible tabs and buttons will be grayed out.

### Adjusting Display Layout – GUI Tab

AlterG has set the default layout in the factory to the first setting listed. Alternative layouts can be made to the navigation bar, Session Metrics, and the Treadmill controls. Below is an example of a selected session layout.

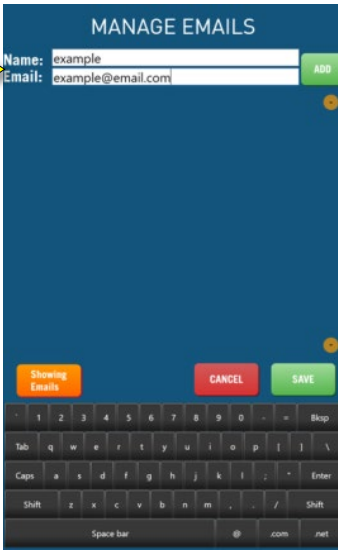
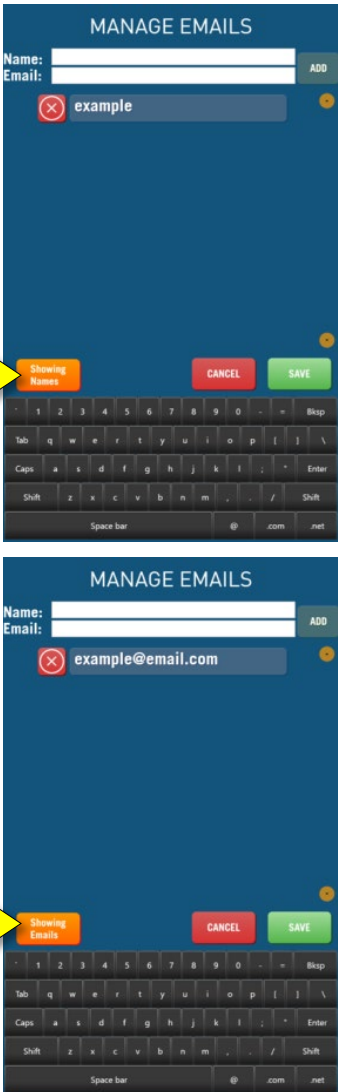


To select an alternative layout, tap the **Change Layout** button adjacent to the layout you want.

You can also change from Metric or Standard units by toggling the **Use Metric Units** button.

Manage Emails – GUI Tab


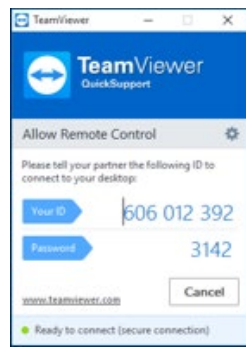
Operator or User Action	Screen
<p>Shorten the End-of-Session process by pre-programming your staff’s email addresses. This allows you to select email addresses from a pre-programmed list and reduce the amount of time spent entering data.</p> <ol style="list-style-type: none"><li>1. Select the GUI tab.</li><li>2. Tap <b>Manage Emails</b>.</li></ol>	A screenshot of a software interface with a dark blue background. At the top, there is a navigation bar with tabs: "GUI", "Support", "Service", "Video", "Registry", and "Diag". Below the navigation bar, there are several buttons: "Manage Emails" (green), "Export Session Data" (white with a red "OFF" toggle), "Use Metric Units" (white with a red "OFF" toggle), "Change Layout" (blue), "Select Logs" (green), "Remove Logs" (green), and a "Return" button (green) in the bottom right corner. A yellow arrow points to the "Manage Emails" button.

Operator or User Action	Screen
<p>3. Enter the clinician's name and email address and tap <b>ADD</b>.</p>	
<p>4. You have the option of displaying the pre-programmed list by name or by email.</p>	

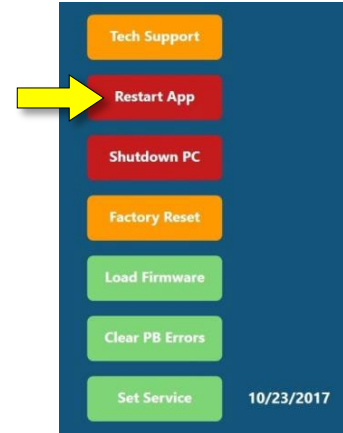
## Tech Support – Support Tab

If you are experiencing issues with your Stride Smart system and are connected to Wi-Fi, our Service Team may be able to get remote access and help you troubleshoot your unit.

Please call our Service Team at (510) 270-5369 if you would like to grant them remote access. They can walk you through the procedure below.


Operator or User Action	Screen
1. In Facility Settings under the Support Tab, tap <b>Tech Support</b> .	
Within 5 seconds you should see a TeamViewer window pop up. After another 5-10 seconds, the ID and password fields will auto-populate.	
2. Read the ID and password back to your technician. This will allow them to connect to your unit and control the screen.	

## Restart App

Operator or User Action	Screen
If you are experiencing issues with Stride Smart, we may ask you to restart your device. To do this, tap <b>Restart App</b> .	


## Shutdown PC

We have designed the Anti-Gravity Treadmill to remain on at all times. The computer and the displays draw exactly 45W, which is less than a 60W light bulb.

Operator or User Action	Screen
<p>Shutdown PC should only be used if Stride Smart and the Anti-Gravity Treadmill must be powered down. Please shut down the unit before unplugging anything from the wall.</p> <p>The system will start up automatically after it is plugged back into the wall.</p>	 <p>The screenshot shows a service menu with several options: Tech Support (orange), Restart App (red), Shutdown PC (red, highlighted with a yellow arrow), Factory Reset (orange), Load Firmware (green), Clear PB Errors (green), and Set Service (green). The date 10/23/2017 is displayed at the bottom right.</p>

## Calibrate Deck Weight – Service Tab

If you are experiencing issues with calibration, you may need to recalibrate the deck weight. Please contact AlterG Service prior to resetting the deck weight.

Operator or User Action	Screen
<ol style="list-style-type: none"><li>1. Select the Service Tab.</li><li>2. To recalibrate deck weight, make sure there is no weight on the treadmill deck. There cannot be anyone standing inside or on the sides of the treadmill.</li><li>3. Tap <b>Set Deck Weight</b>. This number should be around 200lbs (+/- 70).</li></ol> <p>If the calibration issues persist, please contact our AlterG Service Team.</p>	 <p>The screenshot shows the service tab interface. At the top, it says 'COM Port' followed by a dropdown menu showing 'COM11' and a refresh icon. Below this, there are three green buttons: 'Set Deck Weight', 'Pressure Diag', and 'Sim File'. The 'Set Deck Weight' button is highlighted, and the number '230' is displayed next to it. A USB icon is also visible in the top right corner.</p>

## Calibration

You may notice differences in the calibration from older model Anti-Gravity Treadmills. We have implemented a smarter calibration process, which can cause:

- Faster calibration process
- Stronger pressure in the calibration process
- More accurate unweighting
- An inability to unweigh every patient to 20%

The inability to reduce a user's body weight to 20% is the result of increased accuracy during the initial calibration process. If the Anti-Gravity Treadmill is unable to reduce the user's body weight to 20% during calibration, it will only allow them to unweight themselves to the *lowest* accurate body weight point. (that is, if a patient is limited at 37% body weight, this means that during calibration the Anti-Gravity Treadmill could not go below 37%.) A few factors that can cause this, and it should not be of concern unless persistent with all users.

The following may cause an inability to reduce body weight to 20%:

- User moving during the calibration process
- Body composition
- Wearing old shorts (with leaks, tears, or holes)
- Old bag (with leaks, tears or holes)
- Wearing shorts that are not properly fitted to the individual
- Incomplete zipping of shorts to the bag



## Wi-Fi Connectivity

The Anti-Gravity Treadmill is Wi-Fi enabled and can be set up upon installation if the proper Wi-Fi network information is provided.

The Wi-Fi connection can be set up after installation if the facility has a USB keyboard and the Wi-Fi network information. Please contact the AlterG Service Team if you need assistance with connecting to your Wi-Fi network.

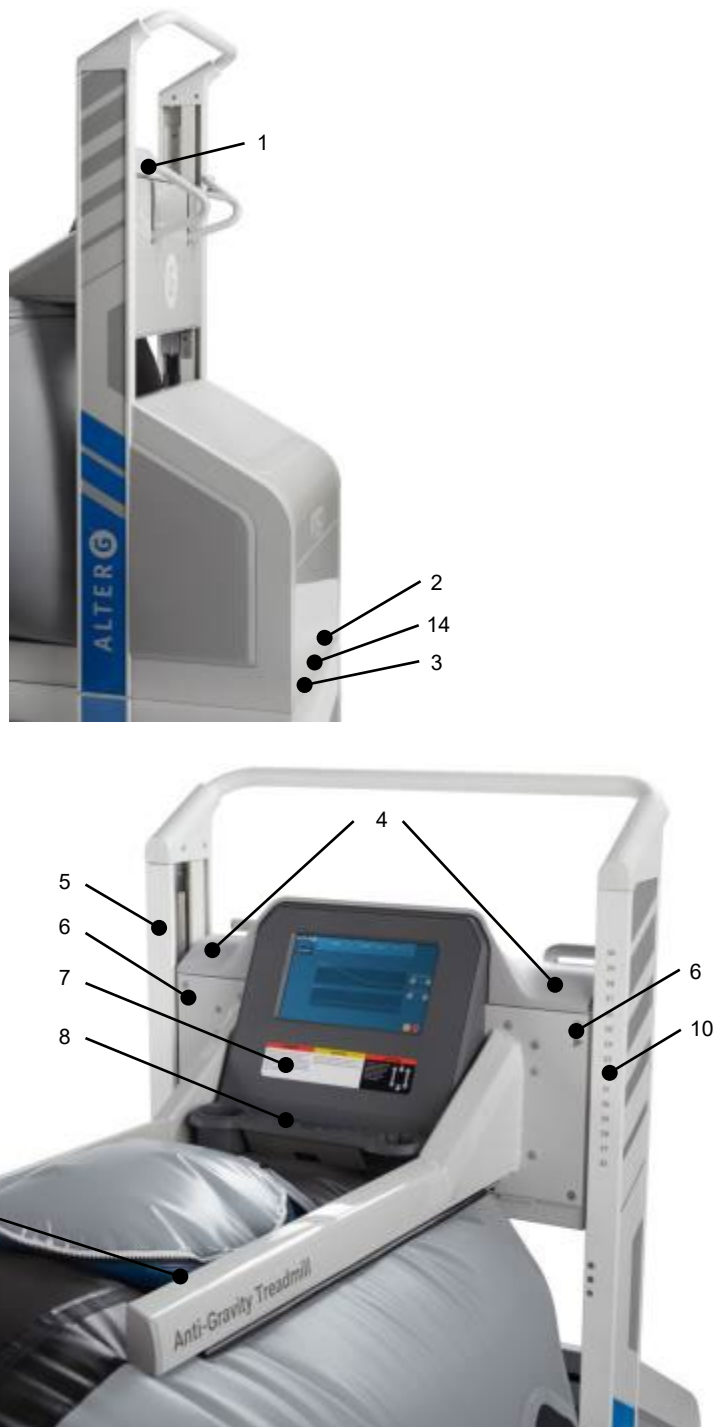
Please note that if your facility has special permissions and security required for Wi-Fi connectivity, we are happy to work with internal IT departments to connect the Anti-Gravity Treadmill to these networks.

## Labels, Locations, Interpretation

You must read and understand the labels on the VIA Anti-Gravity Treadmill. The labels provide information on the operation of the Anti-Gravity Treadmill. Follow all instructions on the labels for a safe and enjoyable exercise experience.

Should any label become damaged and unreadable, contact AlterG immediately to order a replacement.

The locations of the labels are indicated in the diagram below. A graphical representation and detailed description of each label follow.

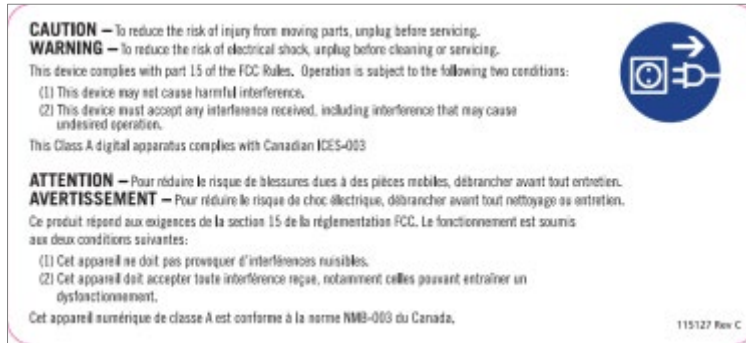


### Label #1



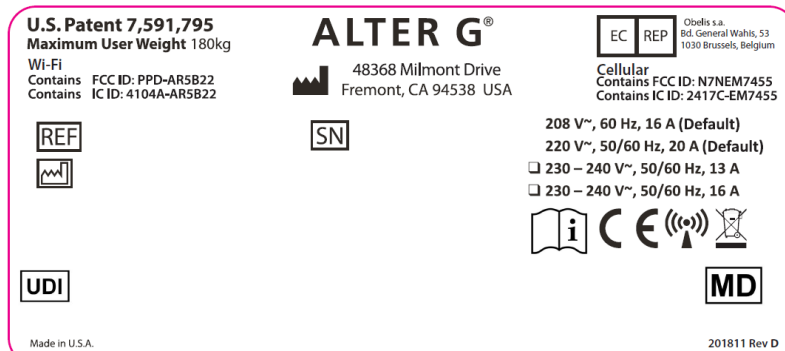
This is the service contact label. It is located on the front of the console above the lockbar and displays the phone number, email address, and website to contact for support and repair.

### Label #2



This label is located on the front panel near the base of the Anti-Gravity Treadmill. Always unplug the treadmill before cleaning or servicing to reduce the risk of injury from moving parts or electrical shock.

### Label #3



AlterG manufactures the system. This label is located on the front panel near the base of the Anti-Gravity Treadmill and identifies the serial and model number of the system, as well as the power and voltage requirements.

### Label #4



This label is located on the display hood, one to the left and one to the right of the touchscreen. Drinks are not allowed on the Anti-Gravity Treadmill anywhere other than in the cup holders. Sealed containers or containers with lids are recommended.

---

### Labels #5 and #10

CM	IN
91	36
	35
86	34
	33
81	32
	31
76	30
	29
71	28
	27
66	26

These labels are located on the vertical stanchions (label #5 indicates centimeters on the left stanchion, label #10 indicates inches on the right stanchion). When a user sets and locks the cockpit, the cockpit height is noted and can be used for future reference.

---

---

### Label #6



This label is located on both sides of the cockpit next to the vertical stanchions and points to the cockpit height set by a user.

---

---

### Label #7



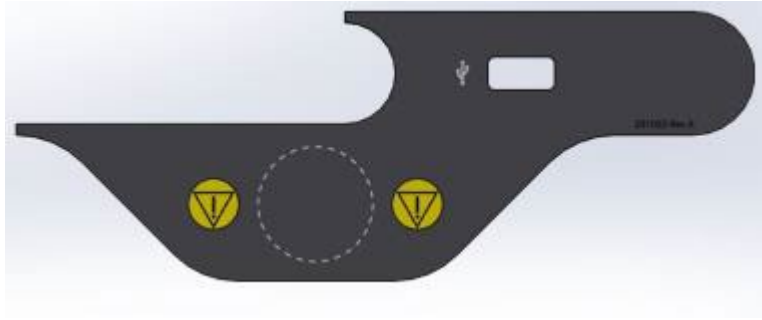
This label is located on the lower right corner of the tubular frame surrounding the touchscreen and console.

The cockpit must be secured in the vertical height adjusters and locked before the user begins exercising. This label is affixed at the front of the cockpit on the top. Be certain to follow these instructions.

---

---

### Label #8



The emergency stop label is located on the control pod below the touchscreen. It indicates where the safety magnet should be placed for operation of the Anti-Gravity Treadmill. If the safety magnet is not placed on the circular stud below the words “EMERGENCY STOP”, the Anti-Gravity Treadmill will not operate. In use, if any sort of emergency should arise, a tug on the lanyard attached to the magnet will displace the magnet and stop the treadmill. Always secure the safety lanyard to the user’s wrist prior to starting the workout.

NOTE: The USB port is intended to be used to transfer End-of-Session data and customized workout programs only. The USB port is not suitable as a charging port for portable devices.

---

---

### Label #9



This label is located inside the handrail on both sides and marks the location of a magnet.

---

---

### Label #11



This label is located within the structure of the Anti-Gravity Treadmill and indicates that a high voltage is present in that location. If you see this label, do not get close to or disassemble any of the components to which it is attached. Because the high voltage can cause serious injury or death, only a qualified AlterG service technician should attempt any repairs.

---

---

### Label #12



This label is located on the inside of the front panel and indicates a connection to earth. Used for zero potential reference and electrical shock protection.

---

### Label #14



**WARNING:** This product can expose you to chemicals including Titanium Dioxide and Carbon Black which are known to the State of California to cause cancer or birth defects or other reproductive harm. For more information, visit [www.P65Warnings.ca.gov](http://www.P65Warnings.ca.gov)

This label is located next to the product label on the front panel near the base of the Anti-Gravity Treadmill. It indicates that this product can expose you to chemicals which are known to the State of California to cause cancer or birth defects or other reproductive harm. For more information, visit [www.P65Warnings.ca.gov](http://www.P65Warnings.ca.gov).

---

## Maintenance

To ensure the safe operation and longevity of your VIA Anti-Gravity Treadmill, you must perform periodic maintenance. You can perform many maintenance tasks yourself; however, it is recommended that an AlterG technician inspect the system every 12 months.



**CAUTION:** Make sure the Anti-Gravity Treadmill is turned off and unplugged before performing any of the maintenance detailed below.

## Disinfection

### Shorts Cleaning and Disinfection

- Before the user puts on or takes off the AlterG's Shorts, make sure they always remove their shoes. Keeping shoes on while putting on the shorts puts a great deal of stress on their seams and will significantly reduce the life of the shorts.
- If any part of the shorts wears out, discontinue use.

**Washing Instructions:** Wash by hand or machine wash on gentle cycle. When using a washing machine place shorts in a mesh bag. Use a mild detergent. Air dry. Do not place the shorts in the dryer.

**Note:** The shorts are latex-free.

Contents:

90% Neoprene

10% Urethane Coated Nylon

AlterG's Shorts should be cleaned and disinfected in accordance with standing clinical policy regarding patient apparel and the degree of exposure risk. Consult the CDC website for the latest guidelines on decontamination of patient equipment and apparel.

Standard cleaning can be carried out with submersion in anti-microbial compounds and mechanical agitation. Follow guidelines for the particular anti-microbial cleaners that are in use when you determine decontamination exposure time and method. AlterG Shorts material construction is of neoprene and urethane; you can consult the manufacturer of your preferred cleaning agents regarding suitability and directions for use. The shorts will tolerate exposure to a 10% bleach solution.

Follow CDC recommended procedures for decontamination when shorts become exposed to human waste or blood or when high-risk patients or high-risk microbial contamination is involved. Disposal of the shorts following exposure to waste, blood or highly contagious microorganisms or when patients at high risk for infection are involved is recommended.

Direct exposure of the shorts to solid waste (feces) blood or broken skin is considered an unusual condition and it may be impossible to adequately disinfect shorts under these circumstances. Shorts which are exposed to higher contamination risk situations should be removed from use beyond the immediate user and sterilized between uses if they are deemed safe for reuse in a particular individual. If broken skin, incontinence or high-risk microbial contamination is possible the situation should be evaluated on a case by case basis.

It is recommended that patients at high risk for urinary or fecal incontinence wear liners, diapers, and other effective means of damming, retention and absorption. Avoid transferring waste or infectious organic matter to the interior of the treadmill because it is very difficult to eliminate organic contamination after it is introduced into the interior of the treadmill.

Urinary catheters and other conduit and bag-based waste storage devices should be used with caution and awareness that the treadmill's internal environment reaches a pressure higher than atmospheric. This pressure difference can pressurize catheter systems, causing them to swell, leak or burst.

## Bag and Frame Cleaning and Disinfection

Wipe the surfaces of the Anti-Gravity Treadmill fabric shell and tubular framework with 10% bleach solution or other detergents/disinfectants that are compatible with urethane coatings and epoxy-based paint films and meet the CDC's guidelines for disinfection. Do not soak surfaces to the point that the bleach solution begins to run. Regular cleaning and wiping of the surfaces after each use is recommended. Following exposure to infectious agents, clean the surfaces of the treadmill in accordance with CDC guidelines or consult AlterG, Inc. **Note:** The bag fabric does not contain latex.

## General Cleaning and Inspection

Periodic cleaning and inspection will help lengthen the life of your Anti-Gravity Treadmill and keep it looking good. The biggest contributor to the failure of the treadmill will be dirt and debris accumulation inside the treadmill. To prevent this, ensure users always wear clean shoes while they exercise. Because the treadmill is a sealed system, the presence of dirt and debris greatly reduces the longevity of the product.

Keeping the system clean will also make it easier to see any problems that might not otherwise be found until it is too late. Below is a general guideline on cleaning and maintenance intervals. If the Anti-Gravity Treadmill is in a dirty environment or under heavy use, cleaning and inspection intervals should occur more frequently. Do not use abrasive brushes or cleaners; they will mar and scratch the paint and plastic surfaces. Also, do not soak any surface with a liquid because the sensitive electronics can be damaged and introduce an electrical hazard.

### Daily

1. Inspect and remove any loose debris from the interior of the bag.
2. Check for abnormal operation. Ensure that there are no unusual performance characteristics such as:
  - Unusual sounds, such as hissing of air leaks, from the treadmill, air blower, or bag.
  - Unusual sights or odors that appear out of the ordinary.
  - Any operational characteristics that have changed, such as reduced treadmill speed or erratic or low bag pressure. Note that low bag pressure can be caused by a miscalibration. Ensure that you have properly followed the calibration steps before you determine that there is a pressurization problem.

### Weekly

1. Check the overall condition of the treadmill.
2. Inspect the bag for tears or leaks.
3. Wipe down exterior surfaces with a damp cloth. This will help prevent the windows from yellowing.
4. Wipe the bag and monitor when needed using a microfiber cloth to avoid scratches. **Note:** Wipe the monitor when the treadmill is off. If you press on the touchscreen when the treadmill is on, you may accidentally activate a function.
5. Check for loose wires and cables.
6. Vacuum the interior of the bag through the access hole in the top of the bag. You can position the cockpit in the highest position and crawl inside the bag for better access.
7. Vacuum around the base of the treadmill.
8. Check shorts for rips or holes.

### Monthly

1. Remove the bag from the frame and thoroughly vacuum the interior of the treadmill. See the instructions that follow.
2. Vacuum any dust accumulated on the screens located inside the blower intake tubes on either side of the front of the treadmill.



## Seasonally

In autumn and winter the drier climate in many regions of the country will cause static buildup when the treadmill is used. Spray the running surface with a staticide spray to prevent static shock to treadmill users and to prevent interference with the treadmill's electronic systems.

## Bag Windows

The windows on your Anti-Gravity Treadmill are made of Strataglass™. Special care must be taken to ensure they remain clean and clear. IMAR™ Strataglass cleaner is recommended for cleaning the clear windows.

**Note:** If your Anti-Gravity Treadmill is exposed to the sun, it is necessary to order a special window cleaner that contains UV protection. For a list of retailers and distributors in your area, contact AlterG or visit Amazon.com and order the IMAR™ Strataglass protective cleaner. If you have any problems with your Strataglass, contact AlterG immediately.

## How Should I Clean My Touchscreen?

When cleaning a touchscreen, it is important to realize that it is sensitive to chemicals, much as is a pair of glasses with plastic lenses (usually polycarbonate with a glare reduction coating). In fact, the cleaning kit supplied with a pair of such glasses is just the ticket for safely cleaning your touchscreen; it typically includes a micro-fiber cloth and a gentle cleaning solution and usually costs about \$10.

### Cleaning Instructions

- Use a soft, lint-free cloth. We especially recommend the 3M Microfiber Lens Cleaning Cloth for cleaning touch panels without needing a liquid cleaner. The material may be used dry or lightly dampened with a mild cleaner or Ethanol. Be sure the cloth is only slightly damp, not wet.
- Never use acidic or alkaline (neutral pH) cleaners, or organic chemicals such as paint thinner, acetone, toluene, xylene, propyl or isopropyl alcohol, or kerosene. Suitable cleaning products are commercially available and pre-packaged for use; one example of such a product is Klear Screen™ or commercially available off-the-shelf retail brands such as Glass Plus® Glass and Surface Cleaner made by Reckitt-Benckiser. Use of incorrect cleaners can result in visual impairment of the touchscreen and damage functionality.
- Never apply cleaner directly to the touchscreen surface; if you spill the cleaner onto the touchscreen, soak it up immediately with an absorbent cloth.
- When using a cleanser, avoid contact with the edges of the film or glass, and with the flex tail.
- Wipe the surface gently; if there is a directional surface texture, wipe in the same direction as the texture.

**Note:** Most products contain 1-3% Isopropyl Alcohol by volume, which is within acceptable limits for Resistive Touch Panel cleaning use.

**Caution:** Many products contain Ammonia, Phosphates, or Ethylene Glycol, which are **NOT ACCEPTABLE**; check product content label carefully.

## Bag Maintenance

- Check the bag for any leaks and note any abnormal hissing sounds. A small leak through the shorts and zipper is normal. If this or other leaks affect the maximum pressure capabilities of the Anti-Gravity Treadmill, contact AlterG for troubleshooting and support.
- Lubricate the zipper on the bag as needed using a zipper lubricant. Always Keep sharp objects away from the bag.
- To provide maintenance on the inside of the Anti-Gravity Treadmill, it may be necessary to remove the bag from the base. It is recommended that this type of maintenance be performed at the end of the day. When the bag is reattached, the foam that forms an airtight seal between the bag and the base requires time to expand and reseal. The procedure for removing the bag for cleaning is described below.

## Bag Removal Procedure for Cleaning

To provide maintenance on the inside of the VIA Anti-Gravity Treadmill, it may be necessary to remove the bag from the treadmill base. Follow the steps in this section to remove, clean inside, and replace the bag.

### Bag Removal and Treadmill Cleaning

1. Power off and unplug the Anti-Gravity Treadmill.
2. You can suspend the bag for better access by securing the bag to the cockpit, and then raising and locking the cockpit.



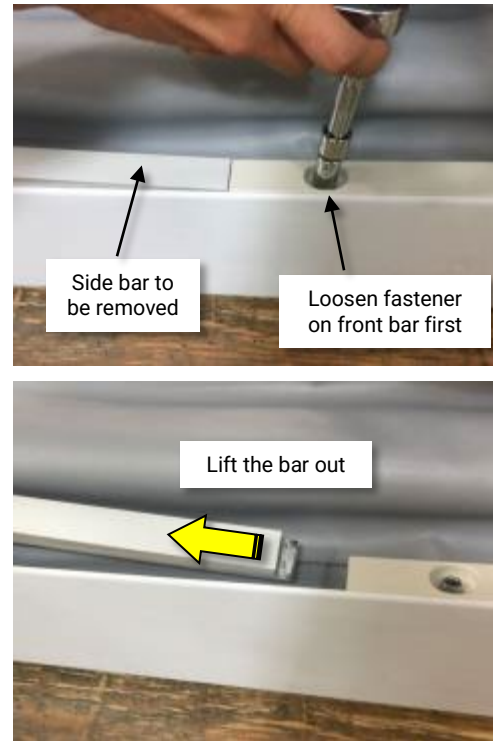
3. Locate the bag-retaining bars on the frame. The bars on the sides and rear are shown. You will remove these three bars. It may be easiest to start with the rear bar.



4. Using a socket wrench with a 10mm socket, unscrew the fasteners from the standoffs.



- 
5. The side bars have a slot on the front end, and the slot fits under the adjacent bar. Loosen the fastener on the adjacent bar, and then lift the side bar out.



- 
6. If a bar is secured tightly, loosen it from underneath using a large screwdriver. Carefully slide the tip of the screwdriver under the bar.



- 
7. Pull the bars completely off the frame and set aside with the fasteners.



- 
8. Lift the base of the bag off the standoffs.

Note that you may not need to remove both the front and back to clean the entire treadmill. Removing the back end and both sides should be sufficient. Typically, it is best to leave the front of the bag attached to the frame.

As shown in the second photo, the bag has been released on both sides and the back. The front of the bag remains attached.



- 
9. Thoroughly vacuum all surfaces you can reach. Get as far under the treadmill as possible.

Clean the back of the treadmill where dirt and debris collect.



- 
10. You can also check the surface of the treadbelt for debris or liquid spills. The treadbelt moves freely as you roll it by hand.

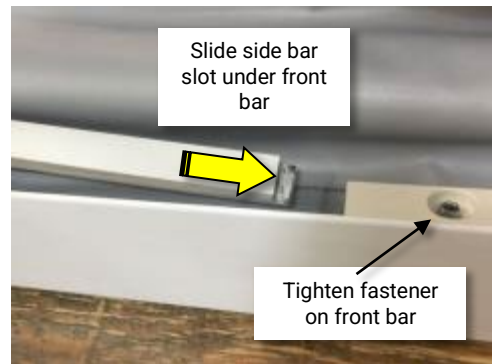


## Bag Reattachment

1. Reattach the bag, lining up the holes with the corresponding standoffs on the frame.



2. Replace the bag-retaining bars on the frame. The side bars must be placed slot-first under the adjacent front bar. Tighten the front bar fastener.



3. Replace and tighten all fasteners. Make sure they are snug. **Do not over-tighten.**



## Appendix A: Anti-Gravity Treadmill Specifications

<b>Performance</b>	<p><b>User Weight Capacity:</b> 80lbs (36kg) – 400lbs (180kg)</p> <p><b>Body Weight Range Adjustment:</b> As low as 20% of user's body weight, in 1% increments (users above 350lbs (159kg) may not achieve 20% unweighting)</p> <p><b>Typical User Height Range:</b> 5'0" (154cm) – 6'3" (190cm)</p> <p><b>Running Surface Area:</b></p> <ul style="list-style-type: none"> <li>• 22" (56cm) wide</li> <li>• 62" (158cm) long</li> </ul> <p><b>Speed Range*:</b> -5 to 15mph (-8 to 24km/hr)</p> <p>*model dependent</p>
<b>Dimensions</b>	<p><b>Length:</b> 85" (217cm)</p> <p><b>Width:</b> 38" (97cm)</p> <p><b>Rubber Mat:</b> 4ft (1.2m) × 9ft (2.7m)</p> <p><b>Weight:</b> 550lbs (249kg), approximately</p>
<b>Recommended Room Dimensions</b>	<p>Provide a footprint of at least 12ft (3.7m) long by 8ft (2.4m) wide for adequate spacing around the treadmill. Leave at least 20 inches (0.5M) of space in front of the unit for access to on/off and electronics panel.</p> <p><b>Note:</b> At the lowest height and high-pressure settings, the bag may expand by as much as ~10" (25cm) per side.</p> <p>Check the ceiling height to ensure that users will not hit their heads on the ceiling while running. The running surface is ~6" (15cm) off the ground.</p>
<b>Environmental</b>	<p><b>Operating Conditions:</b></p> <ul style="list-style-type: none"> <li>• Ambient Temperature : 55°F to 84°F (13°C to 29°C)</li> <li>• Relative Humidity: 20% to 95%</li> </ul> <p><b>Transportation and Storage Conditions:</b></p> <ul style="list-style-type: none"> <li>• Temperature Range: 0°F to 120°F (0°C to 49°C)</li> <li>• Relative Humidity: 20% to 95%</li> </ul>
<b>Electrical Ratings</b>	<p><b>Power Requirements:</b></p> <ul style="list-style-type: none"> <li>• Recommended: 220 VAC 20A, 50/60 Hz</li> <li>• Operational AC Voltage range; 200 - 240 VAC*, 50/60 Hz</li> <li>• Optional AC Voltage range: 120V 20A, 60 Hz*</li> </ul> <p>*At values less than the recommended 220 volts/20A, the ability of the system to reduce body weight to 20% or achieve highest speeds may be compromised.</p> <p><b>Location:</b> Install the front of the treadmill within 12ft (3.7m) from the electrical outlet.</p> <p><b>Electrical Connection (USA):</b> 220V: 20 ampere circuit, NEMA 6-20R receptacle 110V: 20 ampere circuit, NEMA 5-20R receptacle</p> <p><b>International Configuration:</b> The VIA has an IEC 60320 C20 Receptacle for power input. Each unit to come with a mating IEC 60320 C19 Plug and appropriate wall plug for designated country. The appropriate plug should be attached to the power cord of the treadmill using the following wire connection scheme:</p> <ul style="list-style-type: none"> <li>• Blue Conductor: Neutral</li> <li>• Brown Conductor: Line</li> <li>• Ground Conductor: Green/Yellow</li> </ul>



<b>Software utilization, USB access and wireless connectivity</b>	<p>Note: The AlterG Anti-Gravity Treadmill runs on a customized Windows 10 build and is contained within a protected command shell. This prevents a user from exiting the application to access the OS and introduce modifications or hacks. The treadmill can collect Session Data that can be exported from the device. Session Data contains no Personal Health Information (PHI). Session Data can be exported using the USB port where wi-fi is not available. Users should never insert a USB stick of unknown providence into this port. The USB port is compatible with encrypted USB Memory Sticks if added security is required.</p> <p>Session Data is also available for export using a built-in email application (Mailgun) via a wi-fi connection to the product. Mailgun is a secure platform certified to the EU-US Privacy Shield Framework. Mailgun has local servers within both the US and EU and directs traffic to the appropriate server based on the originator's domain.</p> <p>AlterG maintains a MDS2 document available from AlterG regarding secure wi-fi connectivity.</p>
---	--

## Appendix B: Options and Accessories

Please visit the AlterG Store at <http://store.alter-g.com/> or contact your AlterG Sales representative for pricing and ordering.



## Appendix C: Troubleshooting

In most cases, repairs to your Anti-Gravity Treadmill must be completed by an AlterG qualified technician. Contact your AlterG representative, or request repairs at [support@alterg.com](mailto:support@alterg.com).

Before requesting help from a repair technician, you can troubleshoot the problems and potentially resolve them.

### Repairs

Note the following so that we can help you as quickly as possible.

- What is the serial number of the Anti-Gravity Treadmill? The manufacturer's label is located on the treadmill base.
- What happened prior to the problem?
- Did the problem occur unexpectedly?
- Did the problem worsen over time?
- If you hear an unusual noise, from where does the noise originate?
- Was someone using the treadmill at the time the problem occurred?
- Note any other symptoms that might be relevant.
- Does the screen display error messages?

### Air Pressure

If improper pressure is felt during a workout session, check shorts and unit bag for leaks. If pressure issues persist, contact AlterG.

### Treadmill

1. Free wheel: if the treadmill belt is free to move, check and make sure the safety magnet is on the console in the correct location. Next, check that the treadmill is plugged in.
2. The treadmill belt will not move until you start the user interface. If the treadmill belt is moving in another instance, contact AlterG immediately. If the treadmill belt will not move during operation, test other functions, such as incline. If this works, check the screen for any error messages. Record any System Error messages and contact AlterG.

### Leaks

If the fabric bag is torn, or if the shorts are torn, discontinue use and contact AlterG.

### System Errors

The Anti-Gravity Treadmill software has built-in error checking to ensure that all systems are operating within specifications. If an error is detected, "Unexpected Error" is displayed, followed by a description of the detected error. If you see this message, write down the error message and a description of the circumstances under which it occurred.

The error may be the result of an unexpected anomaly that may occur in complex computer-controlled devices. If this is the case, cycle the power from the display console. This may clear the error and correct the problem. If the error persists, contact AlterG. Note the circumstances under which the error occurs and the diagnostic code.

## Appendix D: EMC Statement

**Warning:** The VIA Anti-Gravity Treadmill is MEDICAL ELECTRICAL EQUIPMENT and needs special precautions regarding EMC and needs to be installed and put into service according to the EMC information provided in this manual.

Portable and mobile RF communications equipment can affect the Anti-Gravity Treadmil

The use of accessories, transducers and cables other than those specified by Alter-G Incorporated, may result in increased EMISSIONS or decreased IMMUNITY of the EQUIPMENT.

This EQUIPMENT should not be used adjacent to or stacked with other equipment and that if adjacent or stacked use is necessary, the EQUIPMENT should be observed to verify normal operation in the configuration in which it will be used.

## Appendix E: Warranty

Your Anti-Gravity Treadmill is covered by the following warranty:

- One year parts and labor for the entire machine.

Warranty: AlterG warrants to Customer that the Anti-Gravity Treadmill is free from manufacturing defects for a period of one (1) year from original date of purchase. The Warranty does not cover damage or equipment failure due to misuse, user or other damage, or failure to comply with environmental, electrical requirements and maintenance as outlined in the Anti-Gravity Treadmill User Manual. Any customer modification, disassembly and moving without AlterG Service oversight, or transfer of ownership of the Anti-Gravity Treadmill voids the Warranty and extended Warranty. We invite you to please contact our Customer Success team to allow us to help facilitate the process for you in any circumstance.

Extended Warranty:

AlterG offers an Extended Warranty on a year by year basis for the VIA Anti-Gravity Treadmill as follows:

If you purchase the Extended Warranty at the time of your purchase, AlterG will provide one (1) free preventative maintenance check and service of the VIA by a qualified technician at the end of the first year of use.

An Extended Warranty may be purchased after the sale and installation of the AlterG. For more information, contact your authorized representative or AlterG.

During the Warranty period or Extended Warranty period, AlterG or its authorized service technician will diagnose and repair your VIA Anti-Gravity Treadmill including parts and labor. The service can range from phone calls and emails to onsite service visits as necessary. If you choose not to purchase an Extended Warranty from AlterG, you will be billed at the then current rates for parts and labor plus any travel and/or shipping needed for any service of the product after the initial one (1) year Warranty expires.

Neither the Warranty nor the Extended Warranty covers lost business opportunity due to your VIA Anti-Gravity Treadmill being out of service, nor do the Warranty or the Extended Warranty cover any damage or equipment failure due to misuse and other user damages. This includes: failure to comply with environmental and electrical requirements, as well as the maintenance upkeep protocols outlined in the VIA Anti-Gravity Treadmill User Manual. Any customer modification of the VIA Anti-Gravity Treadmill voids the Warranty. If you must disassemble the VIA Anti-Gravity Treadmill to move it, doing so without an AlterG qualified technician will void the Warranty as well.